

OSTEOMED[®]

A TRADITION OF QUALITY AND INNOVATION
The Logical Curvilinear Solution




LOGIC[™] Mandibular Distraction System

- Low profile intra-oral device for use in cases where distraction osteogenesis is indicated
- Choice of three devices with fixed curvilinear trajectories that mimic the natural growth curve of the mandible
- Utilizes a patented flexible activation wire that provides access through soft tissue for ease of activation
- The activation wire can be cut during the consolidation period reducing the risk of infection while promoting soft tissue healing


LOGIC™

Mandibular Distraction System

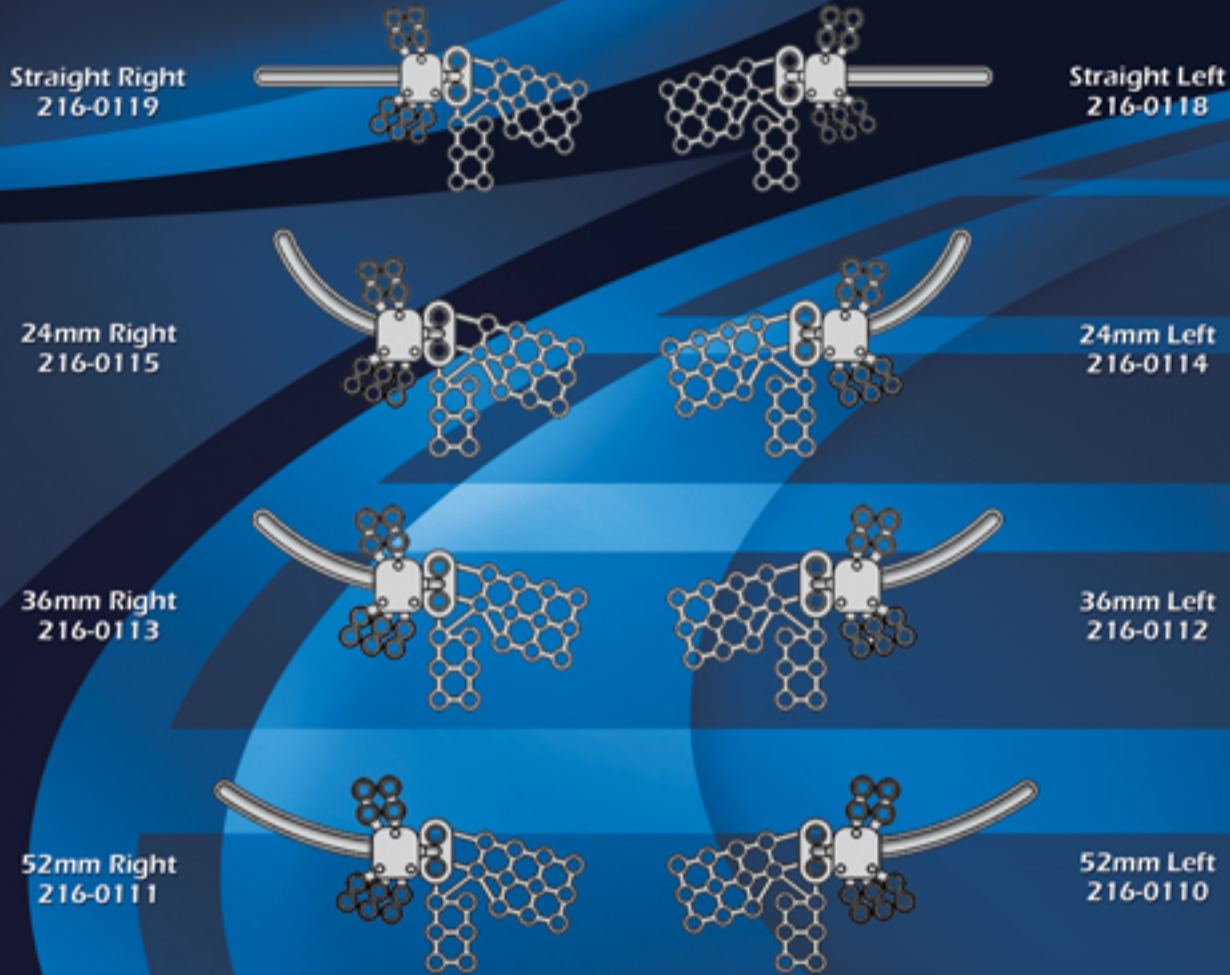

Refer to Product Insert for specific indications and Instructions For Use.

Nickel titanium **activation wires** provide strength and flexibility

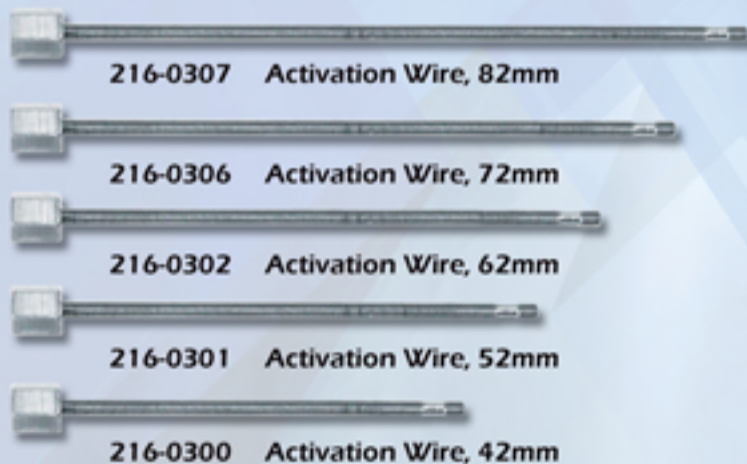

- 216-0307 Activation Wire, 82mm
- 216-0306 Activation Wire, 72mm
- 216-0302 Activation Wire, 62mm
- 216-0301 Activation Wire, 52mm
- 216-0300 Activation Wire, 42mm

2mm and 4mm **spacers** aid in valgus movement



- 216-0227 LOGIC™ Spacer, 2mm
- 216-0228 LOGIC™ Spacer, 4mm



216-0102 Activation Tool

The **activation tool** facilitates rotation of the activation wire.
3 turns = 1mm



A COLSON ASSOCIATE

www.osteomedcorp.com

3885 Arapaho Road • Addison, Texas 75001 • Customer Service: 800.456.7779 • 972.677.4600 • Fax: 800.390.2620

P/N 030-1218 Rev.F